

...better than another set of hands™

Colorectal Kit

Self-retaining retractor system for a wide range of trans-anal surgical procedures


Benefits

- Essential components for key procedures:
 - LAR/Colo-anal Anastomosis
 - Trans-anal tumor excision
 - Anal Fistula Repair
 - Perineal Protectomy (for Rectal Prolapse Repair)
 - Sphincteroplasty
 - Hirschprung's procedures
- Frees up OR personnel to perform other tasks
- Cost effective
- Easy to store, stock, pick and use
- Simplifies surgery
- Provides effective surgical site management with less clutter and better exposure




Lone Star

Colorectal Kit Model # 3707

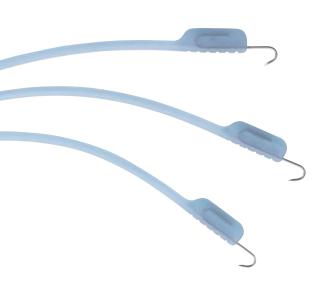



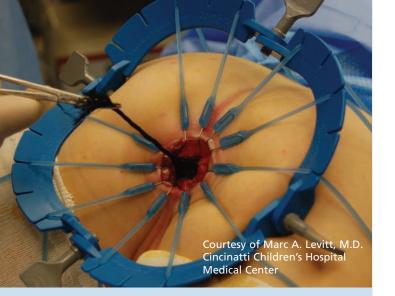

The Lone Star Retractor System facilitates surgical access of the anal canal during Hirschprung's and other colorectal procedures.

Better for the surgeon.
Better for the patient.
Better all around.

Ease and Precision: 3 Simple Steps

1. Place the Retractor Ring around the initial incision. Insert elastic stay hooks in tissue and retract in directly opposing quadrants.

Slots on ring allow easy placement for equal and constant tension.



2. As dissection progresses, reposition hooks and add new ones as needed to achieve full, balanced retraction from any angle.



3. When closing, place stays along axis of the incision to create tension to facilitate precise staple or suture placement.



Lone Star Retractor System™ Colorectal Kit

At last, surgical access the way it should be. Unparalleled exposure. Under your control. This unique retractor system delivers uncluttered access using constant counter-traction with small hooks. The result: a surgical site that you can easily adapt as the procedure progresses.

Get convenience without added cost!


Each kit contains one square disposable retractor ring (14.1 cm x 14.1 cm - Model #3307G) and an eight-pack of 5 mm sharp hook elastic stays (product number 3311-8G) packed in a single-use, sterile dispensing tray.

The cost is the same as if purchased separately, but the convenient packaging facilitates inventory management, making it easy to store, pick and present to the field.

Ordering Information:

PRODUCT NUMBER	DESCRIPTION
3707	Each Lone Star Retractor System Colorectal Kit contains:
	 One square disposable retractor ring, 14.1 cm x 14.1 cm Model# 3307G One eight-pack of 5 mm sharp hook stays – Model# 3311-8G

Each kit is packaged in a single-use sterile dispensing tray. Sold in boxes containing 5 kits.

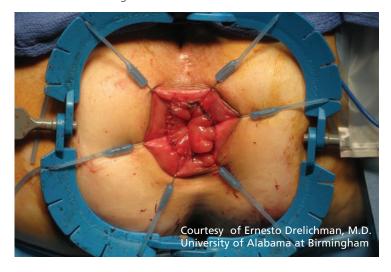

CoperSurgical

CooperSurgical is the leading company dedicated to providing medical devices and procedure solutions that improve health care delivery to women regardless of clinical setting. Our company is fostering that position through expansion of its core businesses and introduction of advanced technology-based products which aid clinicians in the management and treatment of commonly seen conditions.

203.601.5200 | 800.243.2974 | coopersurgical.com | 95 Corporate Dr. | Trumbull, CT 06611

Lone Star Retractor System™ is a trademark of CooperSurgical. U.S. Patent Nos. 5,769,783; 5,785,649; 5,964,697.

© Copyright 2011 CooperSurgical, Inc